



Dolce II Care Toilet Support Arm

Nero

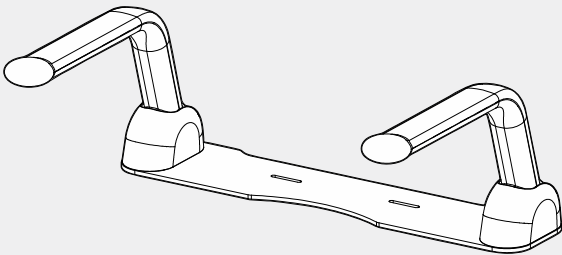

Version : 202607

Drawing

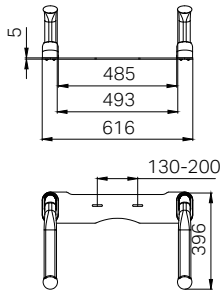

Specifications

SKU
26B002

Warning:

The toilet arms set is designed to be mounted onto toilets. Installation is recommended to be completed by qualified trades person who complies with the AS 1428.1-2021 designed for access and mobility standard. Correct installation is required to ensure safe, long term operation. Incorrect mounting can result in reduced capacity of the rail which can result in injury to user and damage to the wall.

Features


The toilet arms set can be retrofitted onto existing concealed toilets and can be easily removed making them great for rentals and temporary accommodation. Bidet seats can be mounted on top of the arm mounting bracket. 2 sets of fasteners included to suit all types of toilets both concealed and traditional designs. Each arm can fold up out the way when not needed and to give better access from the side of the toilet.

*Dimensions are nominal measurements only.


Installation Instruction

GENERAL INSTRUCTIONS:

- To prevent the lock nut from being pulled out or loosened during use of the arm rest, try to secure the hinge base fixation kit as tightly as possible without causing the nut to be stripped.
- To avoid the Toilet Seat from falling to one side, or being tilted up significantly, do not forcefully shake or put excessive pressure ($C \leq 80\text{kg}$) onto one side of the arm rest after the arm rest is correctly installed in accordance with instructions here.
- To prevent damage to the lifting handrails, the dual handrails simultaneously bear $D \leq 120\text{kg}$.
- The figures below are for reference only.

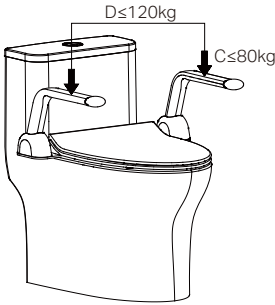


Before Installation:

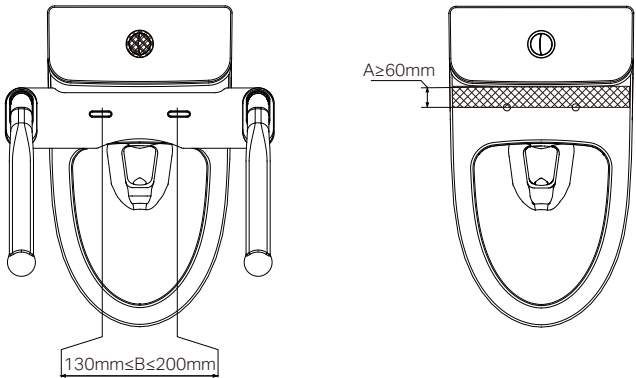


- 1.To prevent interference between the arm rest end and the toilet tank, the area of $A \geq 60\text{mm}$ shall be flat.
- 2.Distance between the installation holes in the toilet seat: $130\text{mm} \leq B \leq 200\text{mm}$.


Top Lock List

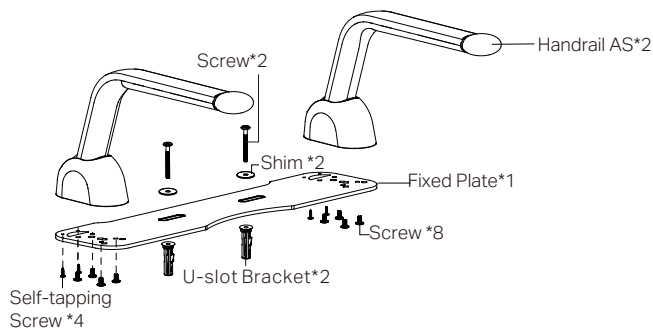

Bottom Lock List

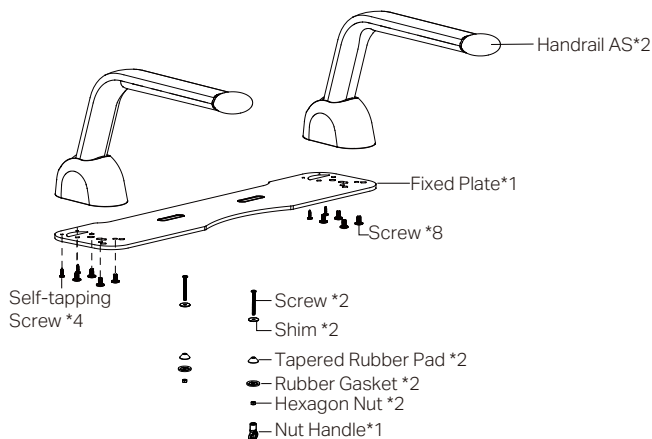


Pre-installation

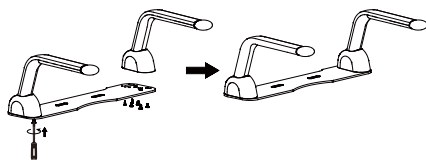


Assemble the handrail. Fix and tighten the screws to assemble both left-hand and right-hand arms to the fixing plate. And depending on your own situation, choose proper one of the next two options for the rest of installation.

Top Lock

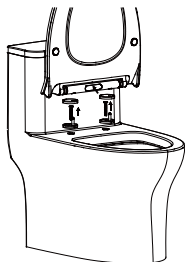


1. Remove the toilet seat and the fixation off the toilet in a proper way.

Note: The above toilet is for reference only. Your toilet and toilet seats may be different.

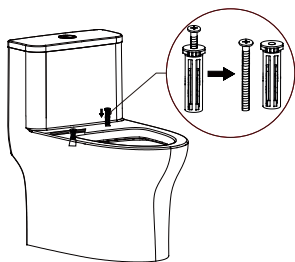


2. Install the handrail screws into the toilet fixing holes.


Note: Do not remove other parts in the integrated screw kit except screws

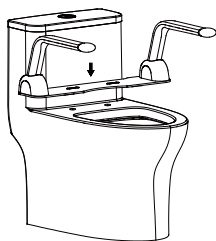


3. Place the fixing plate onto the toilet, align the two installation holes to toilet fixing holes.

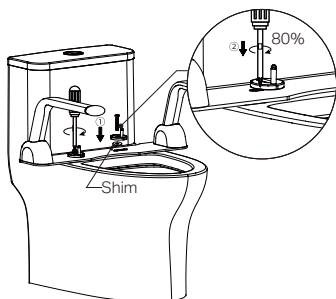

4. Install the fixing plate: install and half-tighten the hinges according to the figure above. (Do not fully tighten the screws)

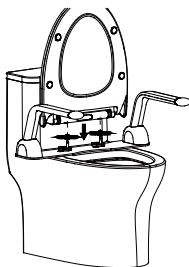


5. Install the toilet seat back.

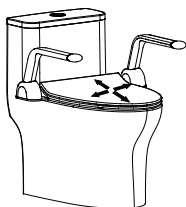


6. Align the toilet seat to the toilet and fully tighten the fixing screws.

(Sometimes it may need to take off the toilet seat to fully tighten the screws).

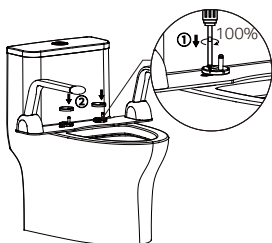

7. Take off the toilet seat to fully tighten the fixing screws, where at the aligning position; Put on the hinge covers.

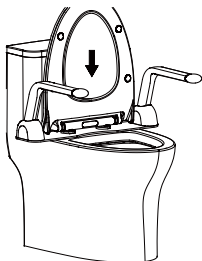


8. Click the toilet seat back.

MAINTENANCE AND CARE:


- Use a mild detergent such as liquid dishwashing soap and warm water for cleaning. Do not use abrasive or high density cleaners that may scratch or dull the surface.
- Carefully read the cleaner product label to ensure the cleaner is safe for use on the material.
- Always test your cleaning solution on an inconspicuous area before applying to the entire surface.
- Wipe surfaces clean and rinse completely with water immediately after cleaner application. Rinse and dry any overspray that lands on nearby surfaces
- Use a soft, dampened sponge or cloth. Never use an abrasive material such as a brush or scouring pad to clean surfaces.

Bottom Lock

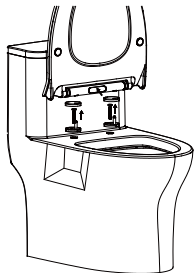


1. Remove the toilet seat and the fixation off the toilet in a proper way.


Note: The above toilet is for reference only. Your toilet and toilet seats may be different.

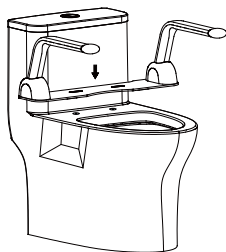


2. Place the fixing plate onto the toilet, align two installation holes to toilet fixing holes.

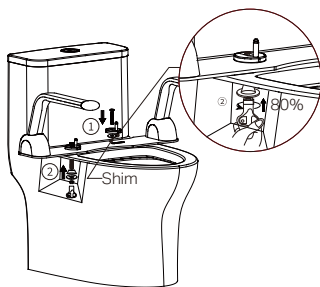

3. Install the fixing plate: install and half-tighten the hinges according to the figure above. (Do not fully tighten the screws).

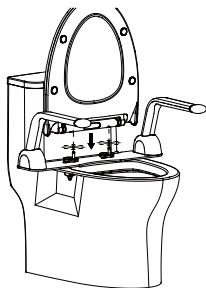


4. Install the toilet seat back.

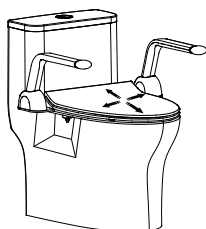


5. Align the toilet seat to the toilet and fully tighten the fixing screws.

(Sometimes it may need to take off the toilet seat to fully tighten the screws).

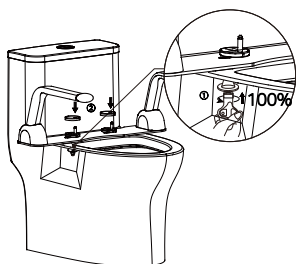

6. Take off the toilet seat to fully tighten the fixing screws, where at the aligning position;
Put on the hinge covers.

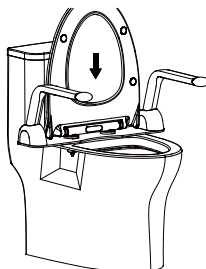

8. Click the toilet seat back.

MAINTENANCE AND CARE:

- Use a mild detergent such as liquid dishwashing soap and warm water for cleaning. Do not use abrasive or high density cleaners that may scratch or dull the surface.
- Carefully read the cleaner product label to ensure the cleaner is safe for use on the material.
- Always test your cleaning solution on an inconspicuous area before applying to the entire surface.
- Wipe surfaces clean and rinse completely with water immediately after cleaner application. Rinse and dry any overspray that lands on nearby surfaces
- Use a soft, dampened sponge or cloth. Never use an abrasive material such as a brush or scouring pad to clean surfaces.

Nero

